

PRODUCT INFORMATION SHEET

Monoclonal antibodies detecting human antigens

CD71

PURE	RUO	REF	IQP-152P	▽	100 tests
FITC	RUO	REF	IQP-152F	▽	100 tests

RUO **For Research Use Only**



Description

Clone DF1513

Isotype IgG1

Specificity Clone DF1513, produces mouse IgG1 immunoglobulins directed against the transferrin receptor. The transferrin receptor is a type II membrane glycoprotein existing as a homodimer of 180-190 kD with the interchain disulphide bond. Clone DF1513 was clustered at the Leukocyte Typing Workshop IV [5].

Antigen distribution

Expression of CD71, or the transferrin receptor, is very low on resting leukocytes but is upregulated upon activation reflecting the iron dependence of activation. In other tissues, CD71 is expressed on most dividing cells and also on brain endothelium and alveolar macrophage [1,2]. The ligand is the serum iron transport protein transferrin.

Summary

CD71 plays a critical role in cell proliferation by controlling the supply of iron, which is essential for many metabolic pathways, through the binding and endocytosis of transferrin, the major iron-carrying protein [1]. The expression of CD71 is regulated at the post-transcriptional level through the control of mRNA stability and is closely linked to intracellular iron levels [4]. Nitric oxide can affect CD71 expression independently of intracellular iron concentration by activating proteins which then stabilize CD71 mRNA [4].

Applications

Monoclonal antibody CD71, clone DF1513 can be applied in flow cytometry, in immunohistochemistry on paraffin sections, cytopots or frozen tissue sections. It is also suitable for use in Western blotting. CD71 antibodies are used as an indicator of proliferation in the typing of lymphomas and leukemias [3].

Usage


All these reagents are effectively formulated for direct immunofluorescent staining of human tissue for flow cytometric analysis using 10 µl/10⁶ leukocytes for singles and 20 µl/10⁶ leukocytes in case of dual and triple combinations. Since applications vary, each investigator should titrate the reagent to obtain optimal results.

HLDA Workshop

Knapp, W., et al eds. 1989 Leukocyte Typing Workshop IV. Oxford University Press.

Representative Data

Clone DF1513 (CD71) was analyzed by flow cytometry using activated peripheral blood lymphocytes from a healthy donor. The cytogram shows direct staining with 10 µl CD71 FITC and 100 µl activated peripheral blood lymphocytes.

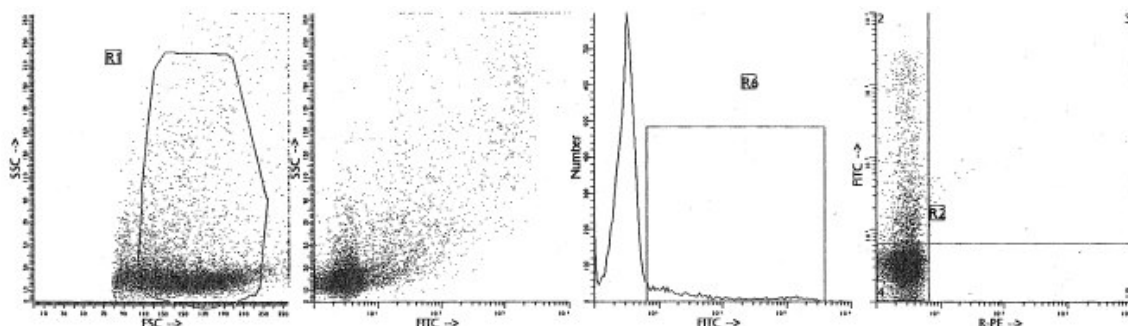

Limitations

1. Conjugates with brighter fluorochromes, like PE and APC, will have a greater separation than those with dyes like FITC and CyQ. When populations overlap, the percentage of positive cells using a selected marker can be affected by the choice of fluorescent label.
2. Use of monoclonal antibodies in patient treatment can interfere with antigen target recognition by this reagent. This should be taken into account when samples are analyzed from patients treated in this fashion. IQ Products has not characterized the effect of the presence of therapeutic antibodies on the performance of this reagent.
3. Reagents can be used in different combinations, therefore laboratories need to become familiar performance characteristics of each antibody in relation with the combined markers in normal and abnormal samples.
4. Reagent data performance is based on EDTA-treated blood. Reagent performance can be affected by the use of other anticoagulants.

Reagents and materials required but not supplied

1. Flow cytometer
2. Flow cytometry disposable 12 x 75-mm capped polystyrene test tubes
3. Micropipette with disposable tips
4. Vortex mixer
5. Centrifuge
6. IQ Lyse - erythrocyte lysing solution (IQP-199)
7. IQ Starfiqs - fixation and permeabilization solution (IQP-200)
8. PBS (phosphate-buffered saline)
9. 1% paraformaldehyde solution in PBS (store at 2-8 °C in amber glass for up to 1 week)

Immunofluorescence staining and lysing protocol

- A - Flow cytometry method for use with purified monoclonal antibodies

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
2. Add to each tube 10 µl of purified monoclonal antibody*. Vortex the tube to ensure thorough mixing of antibody and cells.
3. Incubate the tube for 15 minutes at room temperature in the dark.
4. Wash the labeled cells by adding 2 ml of PBS containing 0.001% (v/v) Heparin, vortexing and centrifuging (2 min 1000 x g.) and discard the supernatant.
5. Add 50 µl of 1:10 dilution of IQ Products F(ab)₂ Rabbit Anti Mouse IgG fluorescent conjugate, [FITC (IQP-190F); R-PE (IQP-190R)] in PBS containing 0.001% (v/v) Heparin to the tube. It is recommended that the tube is protected from light.
6. Mix by vortexing and incubate for 15 minutes at room temperature in the dark.
7. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
8. Incubate for 10 minutes at room temperature in the dark.
9. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
10. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
11. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
12. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

- B - Flow cytometry method for use with labeled (FITC, R-PE, Cy-Q or APC) monoclonal antibodies

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
2. Add to each tube 10 µl of labeled monoclonal antibody*. Vortex the tube to ensure thorough mixing of antibody and cells.
3. Incubate the tube for 15 minutes at room temperature in the dark.
4. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
5. Incubate for 10 minutes at room temperature in the dark.
6. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
7. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
8. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
9. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

- C - Flow cytometry method for use with dual and triple combinations

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
For combinations with anti-kappa and/or anti-lambda Ig see application note below.
2. Add to each tube 20 µl of labeled monoclonal antibody combination*.
3. Vortex the tube to ensure thorough mixing of antibody and cells.
4. Incubate the tube for 15 minutes at room temperature in the dark.
5. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
6. Incubate for 10 minutes at room temperature in the dark.
7. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
8. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
9. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
10. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

* Appropriate mouse Ig isotype control samples should always be included in any labeling study

** PBS: Phosphate Buffered Saline, pH 7.2

Application note for anti-kappa and/or anti-lambda Ig combinations

Add 2 ml of PBS containing 0.001% (v/v) Heparin (**prewarmed to 37 °C**) to the cell suspension
Vortex, centrifuge (2 min at 300x g) and discard the supernatant
Repeat this step twice
Resuspend the pelleted blood cells in 100 µl PBS containing 0.001% (v/v) Heparin



Handling and Storage

Antibodies are supplied either as 100 tests per vial (1 ml) for singles or 50 tests per vial (1 ml) for dual and triple combinations. They are supplied in 0.01 M sodium phosphate, 0.15 M NaCl; pH 7.3, 0.2% BSA, 0.09% sodiumazide (NaN_3). Store the vials at 2-8 °C. Monoclonal antibodies should be protected from prolonged exposure to light. Reagents are stable for the period shown on the vial label when stored properly.

Warranty

Products sold hereunder are warranted only to conform to the quantity and contents stated on the label at the time of delivery to the customer. There are no warranties, expressed or implied, which extend beyond the description on the label of the product. IQ Products is not liable for property damage, personal injury, or economic loss caused by the product.

Characterization

To ensure consistently high-quality reagents, each batch of monoclonal antibody is tested for conformance with characteristics of a standard reagent. Representative flow cytometric data is included in this data sheet.















Warning

All products contain sodiumazide. This chemical is poisonous and hazardous. Handling should be done by trained staff only.

References

1. Testa, U., et al 1993. Crit. Rev. Oncog. 4.241-276
2. Barclay, A.N., et al. 1997. Leucocyte Antigen FactsBook. Academic Press. London
3. Oudemans, et al. 1986 Cancer 58. 1252
4. Hentze, M.W. and Kuhn, L.C. 1996. Proc. Natl Acad. Sci. USA 93. 8175-8182
5. Knapp, W., et al eds. 1989 Leucocyte Typing Workshop IV. Oxford University Press

Explanation of used symbols

	Consult instructions for use
	Catalogue number
	Sufficient for
	In Vitro Diagnostic medical device
	Caution, consult accompanying document
	Keep away from (sun)light
	Biological risks
	Temperature limitation (°C)
	For Research Use Only
	Batch code
	Use by yyyy-mm-dd
	Manufacturer
	Authorized Representative in the European Community
	Conformité Européenne (European Conformity)

		Label - tandem	Ex -max (nm)	Em -max (nm)
P	PURE	purified material	-	-
F	FITC	FITC	488	519
R	R-PE	PE	488, 532	578
C	CyQ	PE-Cy5.18	488, 532	667
A	APC		595, 633, 635, 647	660
PC	PerCP		488, 532	678
PCC	PerCP-Cy5.5		488, 532	695



IQ Products BV
Rozenburglaan 13a
9727 DL Groningen, The Netherlands

 +31 (0)50 57 57 000
 +31 (0)50 57 57 002
 Technical marketing@iqproducts.nl
 Orders orders@iqproducts.nl
 www.iqproducts.nl

 **Products**
bright fluorescence