

PRODUCT INFORMATION SHEET

Monoclonal antibodies detecting human antigens

CD31

R-PE RUO REF IQP-552R ▾ 100 tests | REF IQP-552R50 ▾ 50 tests

RUO

For Research Use Only



Description

Clone MEM-05

Isotype Mouse IgG1

Specificity The antibody MEM-05 reacts with CD31 (PECAM-1), a 130-140 kDa type I transmembrane glycoprotein.

Antigen distribution

CD31 is expressed on monocytes, platelets, granulocytes, endothelial cells and stem cells of the myeloid lineage.

Summary

CD31 is an inhibitory co-receptor involved in rebulation of T-cell and B-cell signaling by a dual immunoreceptor tyrosine-based inhibitory motif (ITIM) that upon associated kinases-mediated phosphorylation provide docking sites for protein-tyrosine phosphatases. CD31 is expressed ubiquitously within the vascular compartment and is located mainly at junctions between adjacent cells. N-terminal Ig-like domain of CD31 is responsible for its homophilic binding, which plays an important role in cell-cell interactions. CD31 is a multifunctional molecule with diverse roles in modulation of integrin-mediated cell adhesion, transendothelial migration, angiogenesis, apoptosis, negative regulation of immunoreceptor signaling, autoimmunity, macrophage phagocytosis, IgE-mediated anaphylaxis and thrombosis. It is one of key regulatory molecules in vascular system.

Applications MEM-05 antibodies are applied in flow cytometry for analysis of blood and bone marrow.

Usage


All these reagents are effectively formulated for direct immunofluorescent staining of human tissue for flow cytometric analysis using 10 $\mu\text{L}/10^6$ leukocytes for singles and 20 $\mu\text{L}/10^6$ leukocytes in case of dual and triple combinations. Since applications vary, each investigator should titrate the reagent to obtain optimal results.

HLDA Workshop

Leukocyte Typing Workshop V

Representative Data

CD31, clone MEM-05, was analysed by flow cytometry using a blood sample from a healthy human volunteer. Direct staining was performed using 10 μl of R-PE conjugated antibody and 100 μl blood sample.

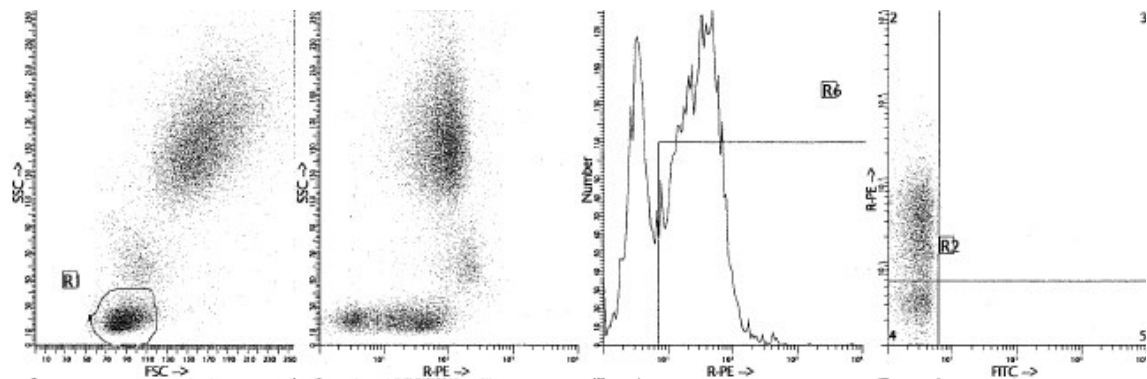

Limitations

1. Conjugates with brighter fluorochromes, like PE and APC, will have a greater separation than those with dyes like FITC and CyQ. When populations overlap, the percentage of positive cells using a selected marker can be affected by the choice of fluorescent label.
2. Use of monoclonal antibodies in patient treatment can interfere with antigen target recognition by this reagent. This should be taken into account when samples are analyzed from patients treated in this fashion. IQ Products has not characterized the effect of the presence of therapeutic antibodies on the performance of this reagent.
3. Reagents can be used in different combinations, therefore laboratories need to become familiar performance characteristics of each antibody in relation with the combined markers in normal and abnormal samples.
4. Reagent data performance is based on EDTA-treated blood. Reagent performance can be affected by the use of other anticoagulants.

Reagents and materials required but not supplied

1. Flow cytometer
2. Flow cytometry disposable 12 x 75-mm capped polystyrene test tubes
3. Micropipette with disposable tips
4. Vortex mixer
5. Centrifuge
6. IQ Lyse - erythrocyte lysing solution (IQP-199)
7. IQ Starfiqs - fixation and permeabilization solution (IQP-200)
8. PBS (phosphate-buffered saline)
9. 1% paraformaldehyde solution in PBS (store at 2-8 °C in amber glass for up to 1 week)

Immunofluorescence staining and lysing protocol

- A - Flow cytometry method for use with purified monoclonal antibodies

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
2. Add to each tube 10 µl of purified monoclonal antibody*. Vortex the tube to ensure thorough mixing of antibody and cells.
3. Incubate the tube for 15 minutes at room temperature in the dark.
4. Wash the labeled cells by adding 2 ml of PBS containing 0.001% (v/v) Heparin, vortexing and centrifuging (2 min 1000 x g.) and discard the supernatant.
5. Add 50 µl of 1:10 dilution of IQ Products F(ab)₂ Rabbit Anti Mouse IgG fluorescent conjugate, [FITC (IQP-190F); R-PE (IQP-190R)] in PBS containing 0.001% (v/v) Heparin to the tube. It is recommended that the tube is protected from light.
6. Mix by vortexing and incubate for 15 minutes at room temperature in the dark.
7. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
8. Incubate for 10 minutes at room temperature in the dark.
9. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
10. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
11. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
12. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

- B - Flow cytometry method for use with labeled (FITC, R-PE, Cy-Q or APC) monoclonal antibodies

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
2. Add to each tube 10 µl of labeled monoclonal antibody*. Vortex the tube to ensure thorough mixing of antibody and cells.
3. Incubate the tube for 15 minutes at room temperature in the dark.
4. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
5. Incubate for 10 minutes at room temperature in the dark.
6. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
7. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
8. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
9. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

- C - Flow cytometry method for use with dual and triple combinations

1. Add 100 µl of EDTA-treated blood (i.e. approx. 10^6 leukocytes) to a 5 ml reagent tube. The content of one tube is sufficient to perform one test.
For combinations with anti-kappa and/or anti-lambda Ig see application note below.
2. Add to each tube 20 µl of labeled monoclonal antibody combination*.
3. Vortex the tube to ensure thorough mixing of antibody and cells.
4. Incubate the tube for 15 minutes at room temperature in the dark.
5. Add 100 µl of IQ Lyse (IQP-199 ready-to-use) and mix immediately.
6. Incubate for 10 minutes at room temperature in the dark.
7. Add 2 ml of demineralized water and incubate for 10 minutes in the dark.
8. Centrifuge the labeled cell suspension for 2 minutes at 1000 x g.
9. Remove the supernatant and resuspend the cells in 200 µl of PBS**.
10. Analyze by flow cytometry within four hours (alternatively, the cells may be fixed by 0.05% of formaline in buffered saline for analysis the next day. Some antigens are readily destroyed upon fixation and this should be taken into account when using this alternative).

* Appropriate mouse Ig isotype control samples should always be included in any labeling study

** PBS: Phosphate Buffered Saline, pH 7.2

Application note for anti-kappa and/or anti-lambda Ig combinations

Add 2 ml of PBS containing 0.001% (v/v) Heparin (**prewarmed to 37 °C**) to the cell suspension
Vortex, centrifuge (2 min at 300x g) and discard the supernatant
Repeat this step twice
Resuspend the pelleted blood cells in 100 µl PBS containing 0.001% (v/v) Heparin



Handling and Storage

Antibodies are supplied either as 100 tests per vial (1 ml) for singles or 50 tests per vial (1 ml) for dual and triple combinations. They are supplied in 0.01 M sodium phosphate, 0.15 M NaCl; pH 7.3, 0.2% BSA, 0.09% sodiumazide (NaN_3). Store the vials at 2-8 °C. Monoclonal antibodies should be protected from prolonged exposure to light. Reagents are stable for the period shown on the vial label when stored properly.

Warranty Products sold hereunder are warranted only to conform to the quantity and contents stated on the label at the time of delivery to the customer. There are no warranties, expressed or implied, which extend beyond the description on the label of the product. IQ Products is not liable for property damage, personal injury, or economic loss caused by the product.

Characterization















To ensure consistently high-quality reagents, each batch of monoclonal antibody is tested for conformance with characteristics of a standard reagent. Representative flow cytometric data is included in this data sheet.

Warning All products contain sodiumazide. This chemical is poisonous and hazardous. Handling should be done by trained staff only.

References

1. Stockinger, H., et al., 1990. J. Immunol., 145 (11): 3889-3897.
2. Simmons, D.L., et al. 1990. J. Exp. Med., 171. 2147-2152
3. Schlossmann et al., eds. 1995. Leucocyte Typing Workshop V. Oxford University Press.

Explanation of used symbols


	Consult instructions for use
	Catalogue number
	Sufficient for
	In Vitro Diagnostic medical device
	Caution, consult accompanying document
	Keep away from (sun)light
	Biological risks
	Temperature limitation (°C)
	For Research Use Only
	Batch code
	Use by yyyy-mm-dd
	Manufacturer
	Authorized Representative in the European Community
	Conformité Européenne (European Conformity)


		Label - tandem	Ex -max (nm)	Em -max (nm)
P	PURE	purified material	-	-
F	FITC	FITC	488	519
R	R-PE	PE	488, 532	578
C	CyQ	PE-Cy5.18	488, 532	667
A	APC		595, 633, 635, 647	660
PC	PerCP		488, 532	678
PCC	PerCP-Cy5.5		488, 532	695



IQ Products BV

Rozenburglaan 13a
9727 DL Groningen, The Netherlands

 +31 (0)50 57 57 000

 +31 (0)50 57 57 002

 Technical marketing@iqproducts.nl

 Orders orders@iqproducts.nl

 www.iqproducts.nl

 **Products**
bright fluorescence